**Antibacterial
Medical Honey™**
**Antibakterieller
Medizinischer Honig™**


Medihoney Pty Ltd
399 Archerfield Rd
RICHLANDS QLD 4077
AUSTRALIA
☎ Toll Free: 1800 006 334
☎ Int Tel: +61 7 3712 8280

Medihoney (Europe) Ltd
200 Brook Drive, Green Park
Reading RG2 6UB UK
☎ UK Toll Free: 0800 071 3912
☎ Int Tel: +44 (0) 118 949 7541

www.medihoney.com

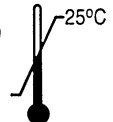STERILE R

109142V1

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Medihoney™ Antibacterial Honey is an exclusive blend of honeys, predominantly *Leptospermum* sp., selected for their unique wound cleaning and antibacterial barrier properties.

Medihoney™ Antibacterial Medical Honey™ is a topical preparation which contains 100% Medihoney™ Antibacterial Honey and can be used for both chronic and acute wound care. It is proven to be clinically effective for:

- providing a moist wound healing environment
- protecting the wound by creating a barrier against wound pathogens including antibiotic resistant strains
- Assisting in achieving a cleaner wound bed
- Rapidly removing malodour
- Reducing the risk of infection from bacteria

Medihoney™ Antibacterial Medical Honey™ is a non-adherent natural dressing and therefore reduces trauma and pain at dressing changes.

Indications

For over-the-counter use, Medihoney™ Antibacterial Medical Honey™ can be used for:

- Superficial wounds such as cuts, scratches, abrasions
- Superficial burns

Under the supervision of a health care professional, Medihoney™ Antibacterial Medical Honey™ can be used for the management of:

- Deep wounds
- Sinus wounds
- Necrotic wounds
- Sloughy wounds
- Infected wounds
- Surgical wounds
- Malodorous wounds

Directions for use

1. Protect wound edges with a barrier cream
2. Medihoney™ Antibacterial Medical Honey™ can be directly applied to the base of the wound. Avoid direct contact between the wound and the nozzle of the tube.
3. Medihoney™ Antibacterial Medical Honey™ can be soaked into a dressing for deep wound cavities.
4. Ensure that Medihoney™ Antibacterial Medical Honey™ is in full contact with the base of the wound – approximately 3mm thickness (1/4 inch). Apply more if wound fluid is present.
5. A non-adherent gauze dressing should be applied to cover the Medihoney™ Antibacterial Medical Honey™
6. A second dressing can be applied if there is a lot of wound fluid. It should be sufficiently absorbent to retain the volume of wound fluid.
7. Medihoney™ Antibacterial Medical Honey™ can be left on the wound for up to 7 days dependent on the amount of fluid present. The dressing should be changed when it is saturated with wound fluid or if the cover dressing is leaking.

Contraindications

- Do not use on individuals who have had an allergic reaction to bee products.

Observations

- When using Medihoney™ Antibacterial Medical Honey™, the wound may initially appear larger in size and depth as unnecessary tissue is cleaned away. This increase should be accompanied by an improved appearance of the wound.
- In deep wounds, packing gauze can be saturated with Medihoney™ Antibacterial Medical Honey™ prior to being introduced into the wound.

Precautions

- For any experience of a reaction other than a short lasting stinging sensation, remove the dressing and wash the affected area.
- For deep, sinus, necrotic, diabetic, sloughy, malodorous or surgical wounds, treatment should be under the supervision of a health care professional.
- Medihoney™ Antibacterial Medical Honey™ may be used on infected wounds under clinical supervision with frequent monitoring of the wound.
- For wounds releasing a high level of fluids, a second dressing of high absorbency is recommended. Not recommended for use on heavily exudating wounds.
- Not recommended for use in body cavities.
- In the event of a slow or non-responding wound please seek medical advice.

Product is only sterile until seal is broken. Do not use if tube is damaged.